

Harpro
seecret



SKIN REJUVENATION

www.seecret.com

AGE IS NO LONGER A SECRET.




COLLAGEN ANTI-AGING

Hapro Seecret® enables your skin to rejuvenate itself continuously in a natural manner. The specially tailored light spectrum stimulates the production of collagen and elastin, which leads to visible improvement of the skin. Discover Hapro Seecret® and experience the sensation of a healthy, attractive appearance. Only you know the secret...


Structure of the skin

Our skin consists of several layers. The protective outer layer is called the Epidermis or scarf skin. Beneath this is a deeper lying layer called the Dermis or true skin. It is this layer that gives the skin its strength; it consists of a dense structure of connective tissue containing collagen and elastin fibres that keep the skin firm and supple. As we grow older, the production of new cells diminishes, reducing the number of collagen and elastin fibres. As a consequence the skin gradually loses its elasticity, resulting in lines and wrinkles.

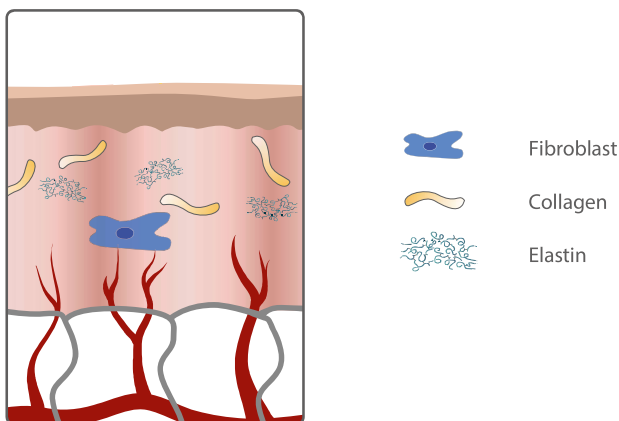

Effect of visible light on the production of collagen

Scientists have observed that exposure of the skin to yellow, orange and/or red light leads to higher production of collagen. This works as follows: In addition to collagen and elastin, the skin also contains fibroblasts. These are special cells that synthesize collagen. Light penetrates the skin to a depth of 10 mm, causing the fibroblasts to grow and divide. Production of collagen and elastin are stimulated in this manner, improving the elasticity of the skin.

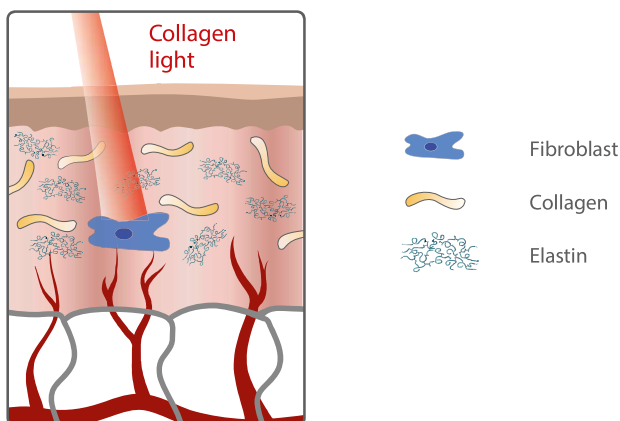


References

Want to know more? The list below contains a number of important scientific publications that prove the effects of collagen light.

- Zhang Y, Song S, Fong C-C, Tsang C-H, Yang Z, Yang M: cDNA microarray analysis of gene expression profiles in human fibroblast cells irradiated with red light. *J. Invest Dermatol.*, 120:849-857, 2003.
- Karu TI, Pyatibrat LV, Kolyakov SF, Afanasyeva NI: Absorption measurements of a cell monolayer relevant to phototherapy: reduction of cytochrome c oxidase under near IR radiation. *J Photochem Photobiol*, 81: 98-106, 2005.
- Lubart R, Friedmann H, Lavi R, Longo L, Jacobi J, Baruchin O, Baruchin AM: A reasonable mechanism for visible light-induced skin rejuvenation. *Lasers Med Sci*, 22:1-3, 2007.
- Lubart R, Lavi R, Friedmann H, Rochkind S.: Photochemistry and photobiology of light absorption by living cells. *Photomed Laser Surg*, 24:179-185, 2006.
- Gao, X., Xing D.: Molecular mechanisms of cell proliferation induced by low power laser irradiation. *J Biomed Sci.*, 16:1-16 online, 2009.

The Light

Approximately 45% of sunlight is visible light we can see with our eyes. In addition to visible light, sunlight also includes components we are hardly able to detect: ultraviolet (UV) and infrared (IR) light. Various scientific investigations have proved that visible light produces a number of responses in the skin. Infrared light, for example, is known to penetrate the skin deeply where, among other things, it improves the blood circulation. Collagen light falls within the colour spectrum of light and is yellow, orange and/or red. The colour spectrum is the bandwidth of light visible to the naked eye. It spans the range of wavelengths from approximately 380 to 760 nm. Collagen light occupies a range of frequencies just below that of infrared light, running from around 600 to 760 nm.

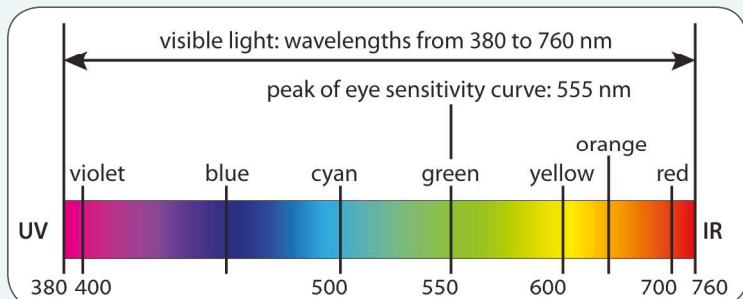


NATURALLY YOUNGER SKIN FROM THE INSIDE OUT

The process of skin regeneration is one of the body's own functions. Hapro Seecret® utilises the effects of pure collagen light, which has a revitalising, stimulating effect. This special light provides extra stimulation to the restorative mechanisms in the skin, resulting in step-by-step, gentle, safe rejuvenation of the skin. Production of collagen and elastin are stimulated, the protein layers are supplemented and the skin's fibre structure is renewed. The improved support and elasticity this provides restores firmness to your skin.

- Hapro Seecret® anti-aging activates the skin's natural production of collagen
- Fine lines and wrinkles fade
- Texture is improved: the skin becomes less oily and feels softer and smoother
- The skin gains a more uniform appearance: age spots fade and your complexion improves

Another positive effect of collagen light is accelerated healing of wounds. Scientific research carried out by the Laser Research Unit, Group of Health Sciences, University of Johannesburg, has demonstrated that light also stimulates cell regeneration, increasing the rate at which wounded skin heals. Moreover, the Hapro Seecret® treatment is a relaxing experience, giving you a pure feeling of physical wellness.

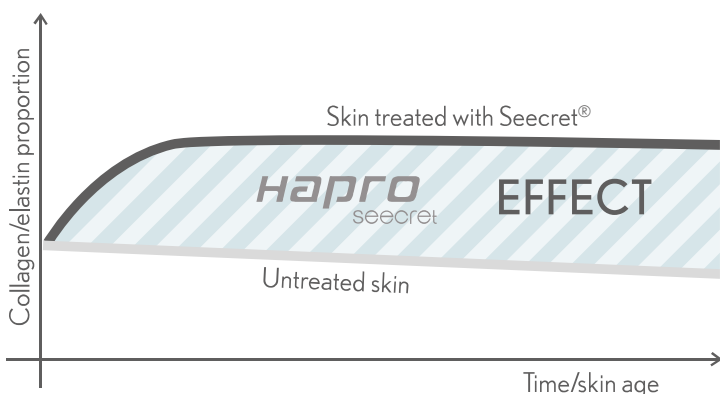

PROF. PAVEL

“The connective tissue of our skin is partly made up of the protein fibres collagen and elastin. One supports our skin, while the other gives it its elasticity. Special cells (called fibroblasts) constantly produce these proteins and thus keep our skin firm, elastic and fresh looking. We are able to regulate this constant reproduction process ourselves up until our mid-twenties. Cell activity decreases as we get older and the skin thus loses its elasticity.”



Prof. Stan Pavel

Associate Professor for Dermatology, former Head of the Phototherapeutic Department and leading researcher of the Dermatology Department of the Leiden University Medical Centre in The Netherlands.

The expert confirms the skin-rejuvenating effect of Hapro Seecret® collagen light treatment.

USAGE

A Hapro Seecret® Session

As with all safe and painless skin stimulation methods, a single Hapro Seecret session will not immediately improve the skin. Aging of the skin is a process that takes place over many years, so it should be no surprise that it cannot be reversed with just one session lasting several minutes. To achieve the desired effects and see visible results a series of rejuvenation sessions is required. We recommend that you indulge in 2-3 rejuvenation sessions per week over a period of 10 to 12 weeks. Keep in mind that a 48-hour rest period is required before you begin the next session.

The results will be measurable after about one month, and you will see visible results after about three months. This differs from person to person, however, and depends greatly on the skin. After three months it may be assumed that the condition of your skin has reached a level at which the frequency of the sessions can be reduced. You can then switch to ‘maintenance sessions’, once every two weeks.

Type	Session time	Number of sessions per week	Number of weeks
Hapro Seecret® C75, C25	30 minutes	2-3	10-12
Hapro Seecret® C700, C300, C200, C180	20 minutes	2-3	10-12

- A course of treatment during a continuous 10 to 12 week period
- 2 to 3 sessions per week
- Rest period of 48 hours between session
- A maintenance session once every two weeks, after completing a course of treatment

As far as is known, from a scientific perspective, there is no need to limit the number of treatments with collagen light. However, based on careful study and practical tests, we established these usage recommendations to help you achieve the most optimal results.

RESULTS

- Reduces fine lines and wrinkles



- Improves skin texture



- Results in more even complexion



UNIQUE AND REFRESHINGLY DIFFERENT

Good, better, best skin improvement – Hapro Seecret® is setting the premium benchmark in this market. Unique details, intelligent technology and a perfect design attract customers and meet all requirements. Hapro Seecret® anti-aging is a comfortable experience with the best results. Both for the professional market as well for use at home. Experience the sensation!



Collagen light

The specially developed Hapro Seecret® collagen light activates the body's own collagen production. The skin will once again produce more collagen and elastin, the two proteins required for a young looking, elastic skin complexion, with fewer wrinkles. Gentle and long-term light treatment with the Hapro Seecret® collagen light will not only rejuvenate your skin in a natural way, but also slow the skin aging process in a long-lasting manner.



Accent light

The additional accent collagen lights, which visibly amplify the anti-aging effect in the head area and on the hands, are unique. The sensitive face and neck areas are noticeably more toned, small lines disappear, deeper contours and wrinkles on the entire body are reduced in a targeted manner.

Welltherm preheating function

Red collagen light is a temperature-neutral light, therefore relatively cold at the start of the session. However, customers are familiar with, and enjoy the comforting warmth of light treatments. As a solution the feature Welltherm provides a warm air flow at the start of the session. It ensures the pleasant feeling of comfort during the Hapro Seecret® anti-aging treatment.



